

TECHNICAL SPECIFICATION

Principles

WBC
Flow Cytometry (FCM) + Laser Light Scatter + Cytochemical Staining Method

RBC/PLT
Electric Impedance Method

HGB
Non-Cyanide Colorimetric Method

Parameters

25 Report Parameters
WBC, Bas#, Bas%, Neu#, Neu%, Eos#, Eos%, Lym#, Lym%, Mon#, Mon%, RBC, HGB, HCT, MCV, MCG, MCHC, RDW-CV, RDW-SD, PLT, MPV, PDW, PCT, P-LCR, P-LCC.

4 Research Parameters
ALY#, ALY%, LIC#, LIC%.

2 Histograms for RBC, PLT
4 Scattergrams for WBC differential

Performance

Parameters	Linearity Range	Precision	Carryover
WBC	0-100x10 ⁹ /L	≤2%(3.5-15x10 ⁹ /L)	≤0.5%
RBC	0-8x10 ¹² /L	≤1.5%(3.5-6.0x10 ¹² /L)	≤0.5%
HGB	0-250g/L	≤1.5%(110-180g/L)	≤0.6%
PLT	0-1000x10 ⁹ /L	≤6.0%(100-149x10 ⁹ /L) ≤4.0%(150-500x10 ⁹ /L)	≤1.0%

Throughput
Up to 60 Samples Per Hour

Display
10.4 Inch Color Touch Display

Data Storage Capacity
Up to 50,000 Results Memory with Scattergrams and Histograms

Communication
LAN Port Supports HI7 Protocol

Interface
4 Pcs USB Ports, 1 Pc Ethernet Port,
Support Uni-Directional LIS
Support Printers, Keyboard, Mouse, Barcode Scanner

Sample Volume
Whole Blood Sample : 20µL
Capillary Whole Blood Mode : 20µL
Pre-Diluted Mode : 20µL

Printout
External Laser Printer/ Inkjet Printer/ Thermal Printer
Various Printout Formats and User-Defined Formats

Operating Environment
Temperature: 10 - 30 °C, Humidity: 70% Max.

Machine Dimension and Weight
Width x Height x Depth : 325mm * 465mm * 470mm
Weight: 30 Kg

Package Dimension and Weight
Length x Width x Height : 630mm * 490mm * 630mm
Weight: 35 Kg

Reagents
Diluent: 20L x 1
Lyse: 60DIFF(500mLx1), 60LH(100mLx1)
Probe Cleanser: 100mL x 2



Copyright © Serachem Diagnostic India Private Limited



SC-Diff 5
AUTOMATED HEMATOLOGY ANALYZER



SC : SERIES OF HEMATOLOGY ANALYZER
Let our experience....
be your guide....

www.serachem-diagnostic.co.in



EC REP
Serachem Diagnostic Limited Serachem Diagnostic, LLC
71-75, Shelton Street, Convent # 354 E, Main Street Suite, 1011 Middletown
Garden, LONDON, WC2H9JQ, D.E. 19709, USA, PH: +1 (302) 601-7242
ENGLAND
E-mail: serachem@serachem-diagnostic.com
www.serachem-diagnostic.com

Manufactured By
Serachem Diagnostic India Private Limited
82, Ward No. 10, Vashisth Colony, Near OHM
Convent School, Kurukshetra, Haryana-136118
E-mail: info@serachem-diagnostic.co.in
www.serachem-diagnostic.co.in
Manufacturing Licence No.: MFG/IVD/2023/000121

Copyright © Serachem Diagnostic India Private Limited

CONTACT US +91 94 1683 1111 E-MAIL US info@serachem-diagnostic.co.in VISIT US www.serachem-diagnostic.co.in

Copyright © Serachem Diagnostic India Private Limited





AUTOMATIC 5 PART HEMATOLOG Y ANALYZER

The Serachem SC-Diff 5 is a compact analyzer that offers fast and high-quality CBC testing for clinics and laboratories. With a 29-parameter complete blood count at the push of a button, the Xpenarray makes on-site hematology easy and affordable.


COMPACT AND POWERFUL

- Space-saving compact design
- Up to 60 samples per hour
- 10.4 inch large touch screen
- 29 parameters + 4 scattergrams + 2 histograms

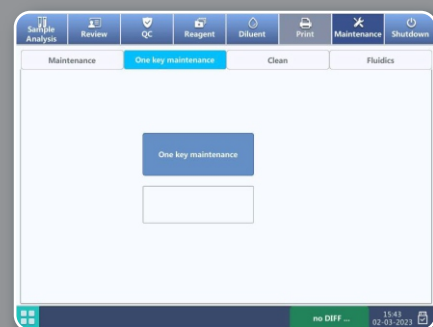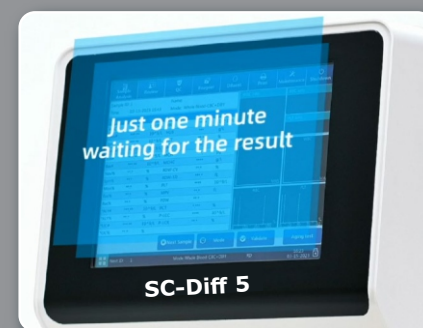


EASY TO USE

- One button analysis
- Cutting-edge technology
- Dual WBC differential channels
- Smart maintenance system

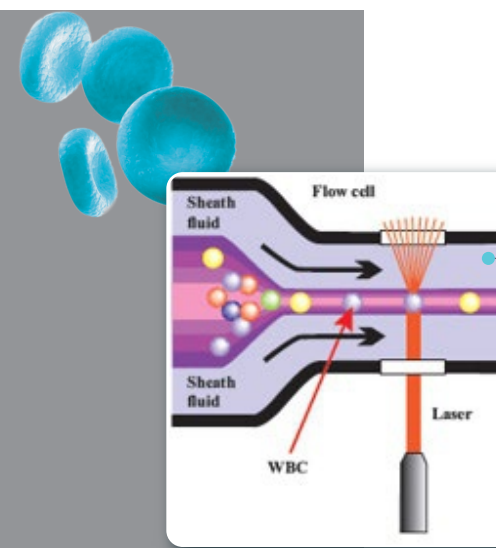


DUAL ACCELERATION FLOW CELL SYSTEM

Reduces the velocity gap between sheath flow and blood cells, providing a more stable signal for accurate measurement via laser scatter

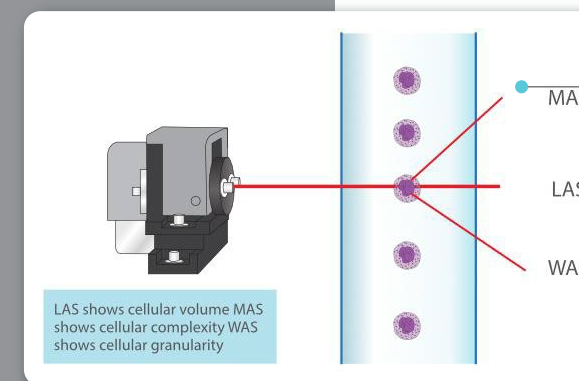

LATEST INNOVATION TRI-ANGLE LASER SCATTERING

The tri-angle laser scattering contributes to more accurate counting. The optical detector receives scatter light signals and converts them into electrical pulses.

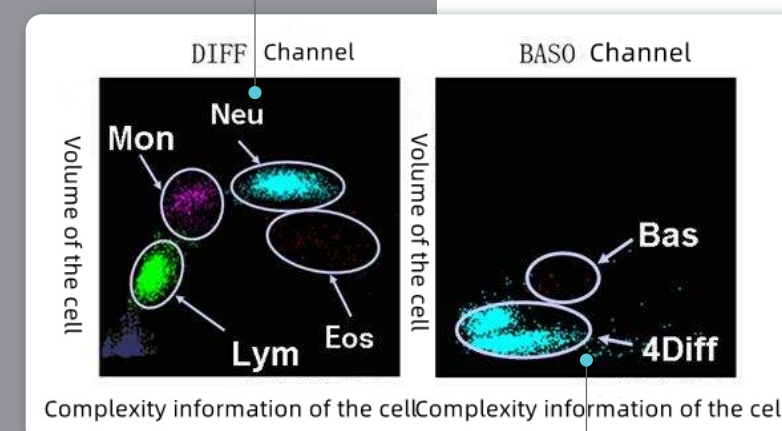

CHEMICAL DYE METHOD BY ORIGINAL REAGENTS

The basophils and eosinophils are differentiated by the unique chemical dye method. The reagents, control and calibrator materials act as one complete system to ensure high accuracy in differentiating between basophils and eosinophils

TWO SEPERATE CHANNELS TO DIFFERENTIATE WBC

Counts lymphocyte, monocyte, neutrophil and eosinophil in "Diff" channel, and basophil in a separate "Baso" channel. Optical method to test Basophil provides more reliable information about cell physiology.